

Otto Bock[®]

QUALITY FOR LIFE



Motus

Product Information





ONE FOR ALL!

Why?

Quite easily, because the adjustment possibilities and options allow to merge with your requirements and needs. And if some of the assemblies seem familiar to you, you are right. We stake on proven technology with new look. Simply, we are offering you our known quality and try to catch your interest.

Do you want to know more ?

MOTUS



The Motus is an adaptive folding wheelchair, provided with a newly designed back upholstery available in two different colors - a real eye catcher. It is an attention getter regardless if you choose a punchy green or muted bronze design.

Crossbrace & Stability



We place great emphasis on the known crossbrace technology. It has proven itself many times and is one of the most successful discoveries in the wheelchair's development. It allows you to manage everyday life situations safely and naturally and fold your Motus for transport in a few easy steps.

Seated in support brackets, the construction offers you a stable driving behaviour in the long term.






Sluggish?

Quite the Contrary!

Motus offers you 48 individual adaptation possibilities by its well thought out rear wheel adapters. It can be configured directly to your body height and driving behavior. Regardless if you prefer active driving or anti-tip safety. You can determine center of gravity which is optimal for you and change it again, if need be. Motus offers you a great variability and adaptability according to your needs!

Your sitting posture can be especially adjusted to your everyday situations and home environment. The same applies to seat height adjustment allowing to drive under the desk, or to seat angle adjustment for ensuring stability as much as possible.

Motus is an all-rounder!

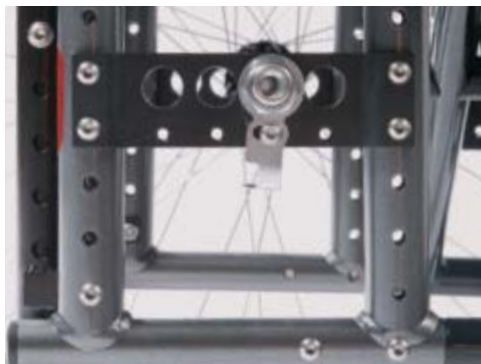

Nice and Practical!

The back padding with an integrated three-dimensional knit spacer fabric is not only a trendy style but it also makes something good for your body. This material allows air to circulate and simultaneously prevents the wind from blowing through. The padding feels pleasantly soft on the skin. The adaptable back upholstery is perfectly adjustable with hook & loop straps and delivered as standard in the basic model as we know that the optimal sitting posture is very important for you!

The padded armrests can be adjusted in height and go very well with Motus. They can be folded away to the rear for easy transfer. This is where quality, design and functionality meet.

The plug-on side panels with armrest pad are also normally available as an option for space-saving transport.



A Pleasing Rest?

How?

You can choose the optional elevating legrests and give rest phase to your legs. The technology allows you to lock the desired position in a few easy steps. Then again, you release the locked height adjustment using only one hand. Easy handling and comfort.

Of course the footrest versions offer angular setting for precise positioning of the feet according to your well-being needs.





In a Nutshell

Frame material	aluminum
Weight	from approx. 13 kg / 29 lbs
Concept	crossbrace folding wheelchair
Center of gravity adjustment	48 positions adjustable through the rear wheel adapter
Brakes	wheel locks, drum brake
Maximum load capacity	125 kg / 275 lbs (seat width up to 48 cm / 19 in) 140 kg / 309 lbs (seat width up to 50.5 cm / 20 in)
Back	straight, height adjustable / 30° ankle adjustment
Back height	30 – 50 cm / 11.9 - 19.7 in
Seat width	35.5 – 55.5 cm / 14 – 21.9 in
Seat depth	36 – 54 cm / 14.2 – 21.3 in
Lower leg length	34 – 52 cm / 13.4 – 20.5 in
Overall width:	52.0 – 72.0 cm / 20.5 – 28.3 in with standard rear wheels * 55.5 – 75.5 cm / 21.9 – 29.7 in with rear wheels with drum brake *
Overall length:	80.5 – 108.5 cm / 31.7 – 42.7 in with 22-in rear wheels 83.0 – 111.0 cm / 32.7 – 43.7 in with 24-in rear wheels
Frame colors	silver metallic, anthracite metallic


* Applicable to push ring attachment, narrow (for push ring attachment, wide: + 20 mm / 0.79 in) and a 0° camber of the rear wheels.



Frame color: anthracite metallic



Frame color: silver metallic

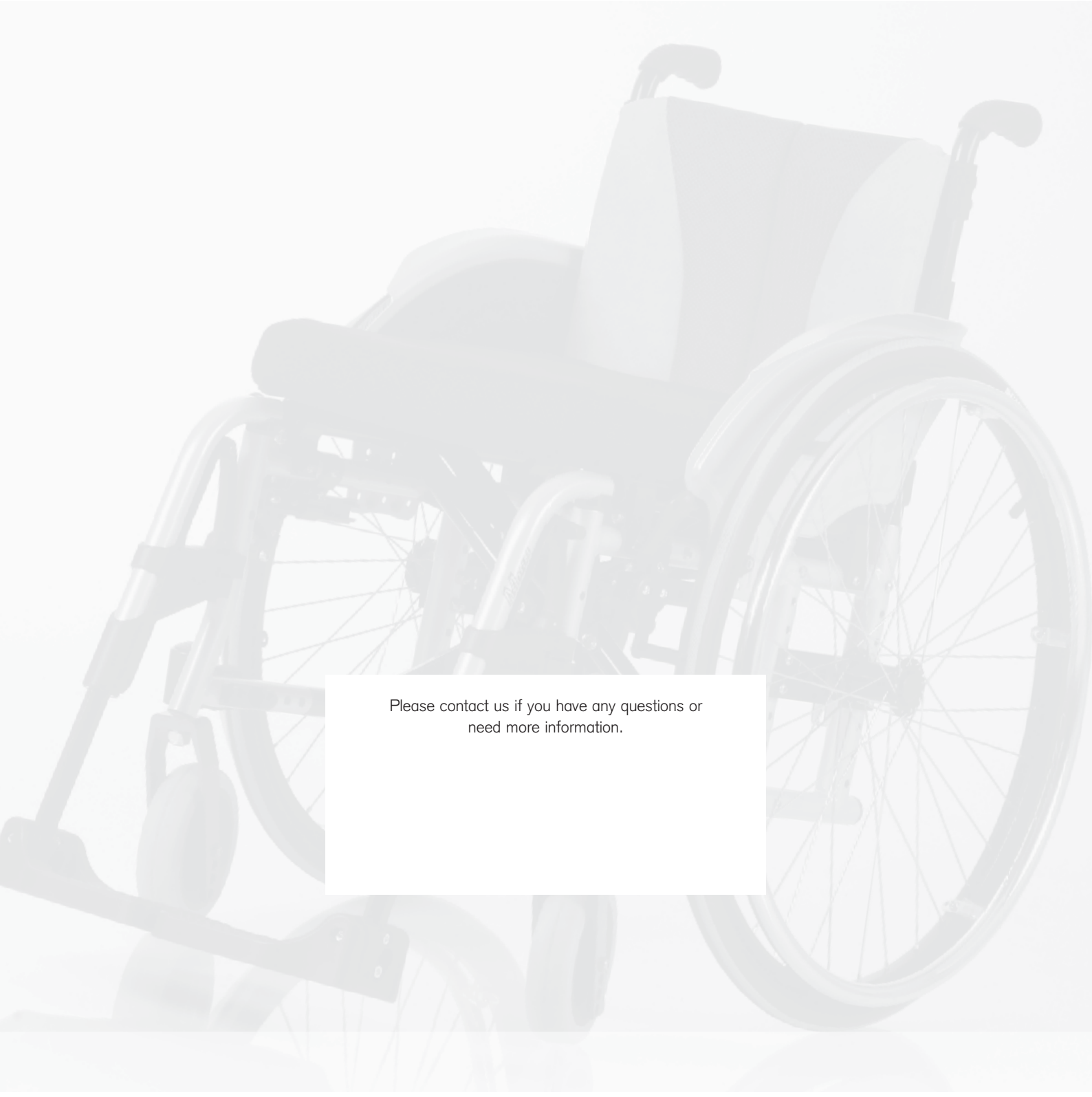

Please contact us if you have any questions or
need more information.

Otto Bock[®]

QUALITY FOR LIFE

Otto Bock HealthCare GmbH

Max-Näder-Straße 15 · 37115 Duderstadt/Germany · Phone +49 5527 848 1590 · Fax +49 5527 848 1676 · reha-export@ottobock.de · www.ottobock.com